


# The SpineJet® Endoscopic Discectomy System

*Experience Fluidjet  
Power First Hand*

*Combining the power of  
fast, precise, controlled  
nucleus removal with  
enhanced visualization*

**Designed to work through any  
endoscope with a 3.1mm  
working channel providing  
predictable nucleus  
removal regardless  
of disc hydration**

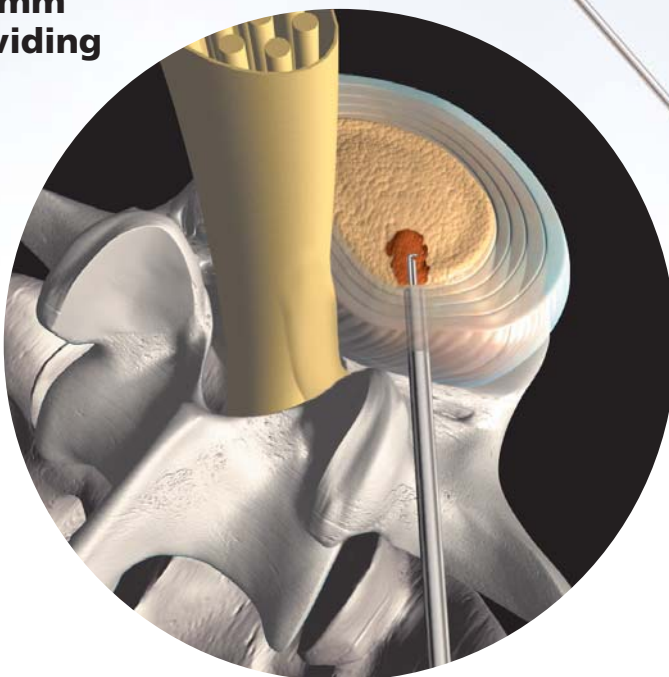

MEDICAL  
DESIGN  
EXCELLENCE  
AWARDS®  
2006 WINNER

HydroCision®  
The Leader in HydroSurgery®

## Efficient


- Simultaneously cuts and aspirates nucleus reducing procedural time
- Power driven controlled nucleus removal
- Sterile, disposable handpiece for faster set-up eliminating the need to clean, inspect, and sterilize before each procedure
- Reduces need for multiple conventional instruments

## Safe

- Cool, fluid based cutting action eliminates risk of thermal tissue damage
- Round atraumatic EndoResector tip design reduces risks of annular puncture and endplate damage
- Instrument insertion depth gaged by 0.5cm laser markings
- Depth indicator controls instrument insertion depth
- Reduces instrument passes through the surgical site minimizing the risk of inadvertent nerve injury

## Effective

- Removes precise, predictable amount of nucleus to decompress the nerve root
- Ablates tissue efficiently regardless of patient's age and disc hydration

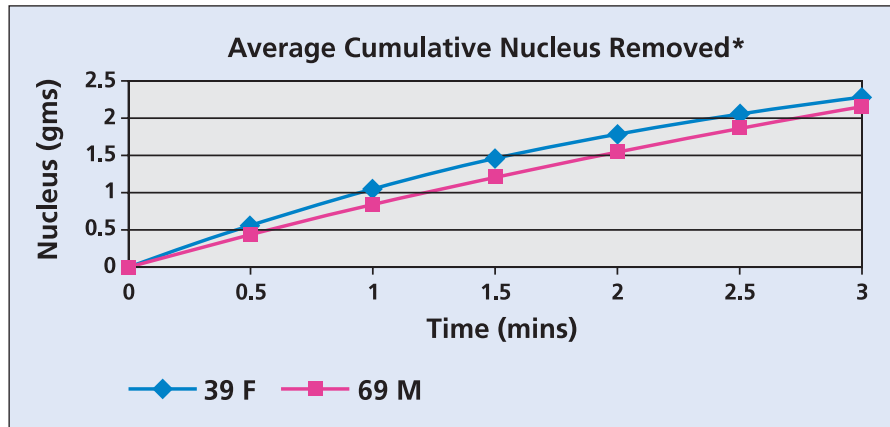

Removes approximately 2.3 grams of nuclear material in a three minute run time regardless of disc hydration.

\* Cadaver studies. Data on file.  
HydroCision, SpineJet, and The Leader in HydroSurgery are registered trademarks of HydroCision, Inc.  
©2004-2006. All rights reserved. Patent No.6,375,635, 6,960,182, 6,923,792, 7,122,017 other US and Foreign Patents Pending.  
P/N 1000-2460 Rev B 07/09 TG-9922

## Ordering Information

SpineJet Endoscopic Discectomy System	
Description	Catalog No.
SpineJet EndoResector Handpiece <i>(disposable, supplied sterile)</i>	56110
SpineJet Quick Connect <i>(Quick Connect includes tubing set - disposable, supplied sterile)</i>	55400
SpineJet QuicKart HydroSurgery System <i>(Includes SpineJet HydroSurgery Console and QuicKart Procedure Cart)</i>	50370
SpineJet HydroSurgery Console <i>(Includes Foot Pedal and Power Cord)</i>	52700
QuicKart Procedure Cart	50300

## HydroCision Customer Care Line 888.747.4470

Call Toll Free for Product Support  
267 Boston Road, Suite 28  
North Billerica, MA 01862-2336  
Tel. 978.474.9300  
Fax 978.600.5058  
Email: [info@hydrocision.com](mailto:info@hydrocision.com)  
[www.hydrocision.com](http://www.hydrocision.com)