


# The SpineJet® Percutaneous Hydrodiscectomy System

*Experience Fluidjet  
Power First Hand*

***Precise, Predictable  
Nucleus Removal***

***Power driven controlled nucleus  
removal for decompression  
regardless of disc  
hydration***

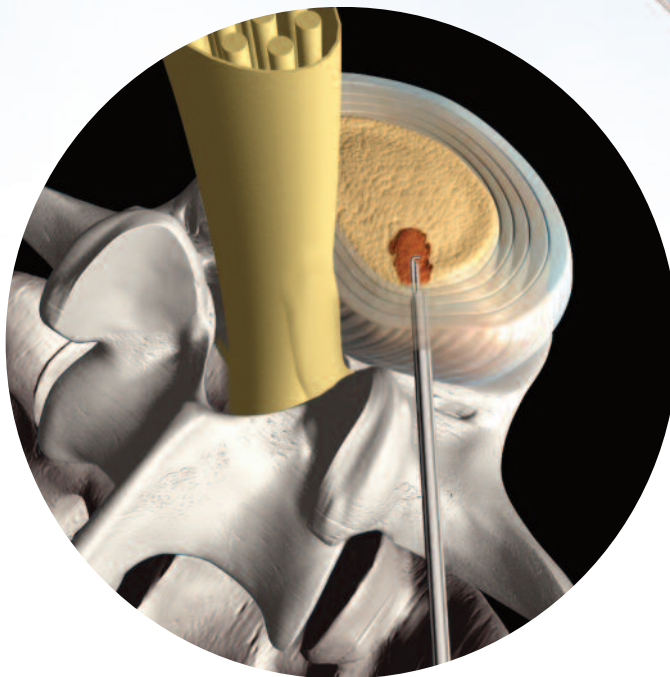

Learn more at [www.WashAwayBackPain.com](http://www.WashAwayBackPain.com)

**HydroCision®**  
*The Leader in HydroSurgery®*

## Efficient

- Simultaneously cuts and aspirates nucleus reducing procedural time
- Power driven controlled nucleus removal
- Sterile, disposable access set and handpieces for faster set-up, eliminating the need to clean, inspect and sterilize before each procedure


## Safe

- Cool, fluid based cutting action eliminates risk of thermal tissue damage
- Round atraumatic PercResector tip design reduces risks of annular puncture and endplate damage
- One step seamless insertion of dilator and cannula ensures atraumatic, smooth entry into the disc space dilating, rather than cutting, annular fibers

## Effective

- Removes precise, predictable amount of nucleus to decompress the nerve root
- Ablates tissue efficiently regardless of patient's age and disc hydration
- Small dilated annulotomy reduces risk of reherniation


The HydroCision Percutaneous Hydrodiscectomy System includes both a straight and a curved access cannula.

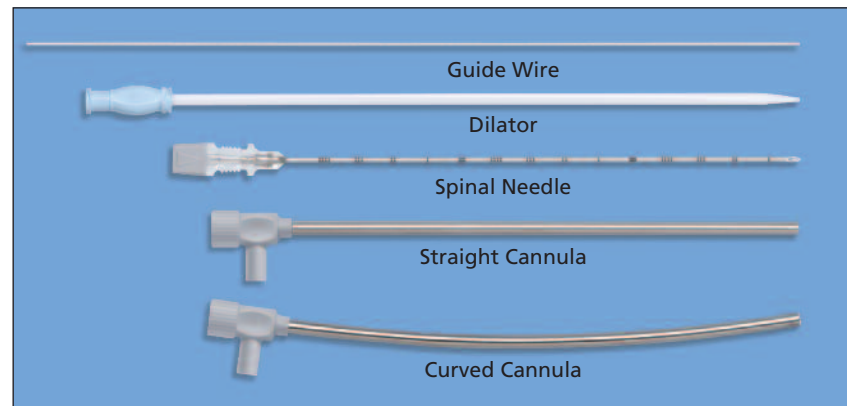

## Ordering Information

### SpineJet Percutaneous HydroDiscectomy System

Description	Catalog No.
<ul style="list-style-type: none"> <li>• SpineJet PercResector Handpiece and Direct Connect Tubing Set</li> <li>• SpineJet Percutaneous Straight and Curved Access Set Options</li> </ul> <i>(Includes one access set which can be configured as straight or curved. Set includes one each straight and curved cannula, guide wire, spinal needle with stylette, and dilator. Disposable, supplied sterile.)</i>	56000
SpineJet HydroSurgery Console <i>(Includes Foot Pedal and Power Cord)</i>	52700
QuickKart Procedure Cart	50300

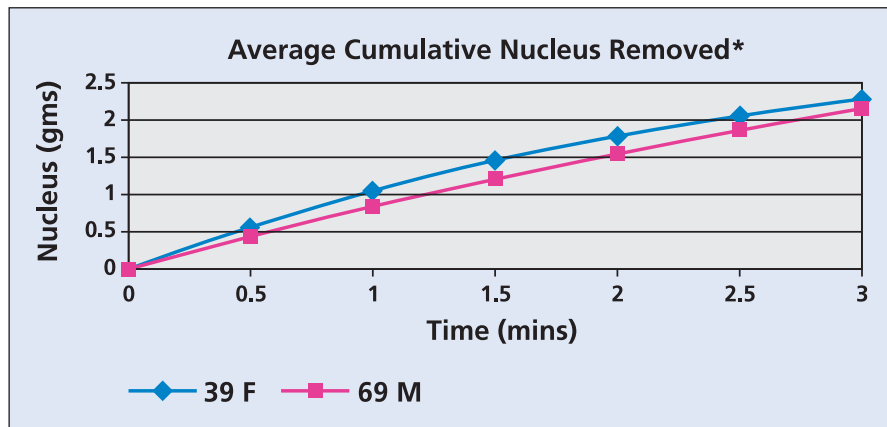

Removes approximately 2.3 grams of nuclear material in a three minute run time regardless of disc hydration.



## HydroCision Customer Care Line 888.747.4470

Call Toll Free for Product Support  
267 Boston Road, Suite 28  
North Billerica, MA 01862-2336  
Tel. 978.474.9300  
Fax 978.600.5058

Email: [info@hydrocision.com](mailto:info@hydrocision.com)  
[www.WashAwayBackPain.com](http://www.WashAwayBackPain.com)

HydroCision, SpineJet, and The Leader in HydroSurgery are registered trademarks of HydroCision, Inc.  
©2004-2006. All rights reserved. Patent No.6,375,635, 6,960,182, 6,923,792, 7,122,017 other US and Foreign Patents Pending.  
P/N 6000-0004 Rev A (replaces 1000-1658) 3/12 KAM189